



Oris™ Cell Migration Assay

Product No.: CMA1.101 & CMA5.101

96-well, 2-D Assay for Investigating
Cell Migration of Adherent Cell Lines

PROTOCOL & INSTRUCTIONS

Table of Contents

I.	INTRODUCTION.....	2
II.	ORIS™ PLATE DIMENSIONS	2
III.	MATERIALS PROVIDED	3
IV.	MATERIALS REQUIRED	3
V.	CELL MIGRATION ASSAY PROTOCOL.....	4
VI.	DATA ACQUISITION	6
VII.	ORDERING INFORMATION.....	8
VIII.	TERMS & CONDITIONS.....	8
	APPENDIX I - Determining Optimal Cell Seeding Concentration.....	9
	APPENDIX II - Determining Optimal Fluorescence Plate Reader Settings...9	



AMSBIO | www.amsbio.com | info@amsbio.com



UK & Rest of the World

184 Park Drive, Milton Park
Abingdon OX14 4SE, U.K.
T: +44 (0) 1235 828 200
F: +44 (0) 1235 820 482



North America

1035 Cambridge Street,
Cambridge, MA 02141.
T: +1 (617) 945-5033 or
T: +1 (800) 987-0985
F: +1 (617) 945-8218



Germany

Bockenheimer Landstr. 17/19
60325 Frankfurt/Main
T: +49 (0) 69 779099
F: +49 (0) 69 13376880



Switzerland

Via Lisano 3
(CP.683)
6900 Massagno
T: +41 (0) 91 604 55 22
F: +41 (0) 91 605 17 85

Oris™ CELL MIGRATION ASSAY


I. INTRODUCTION

The Oris™ Cell Migration Assay is a reproducible, sensitive, and flexible assay that can be used to monitor cell migration. Formatted for a 96-well plate, the assay utilizes Oris™ Cell Seeding Stoppers made from a medical-grade silicone to restrict cell seeding to the outer annular regions of the wells. Removal of the stoppers reveals a 2mm diameter unseeded region in the center of each well, i.e., the detection zone, into which the seeded cells may then migrate. The Oris™ Detection Mask is applied to the plate bottom and restricts visualization to the detection zones, allowing only cells that have migrated to be detected (see Figure 1). The Oris™ Cell Migration Assay is designed to be used with any commercially available stain or labeling technique. Readout can be performed by microscopy or use of a microplate reader.

The Oris™ Cell Migration Assay system has been designed for use with adherent cell cultures. This assay has been successfully used with HT-1080, PC-3, A549, NCI H1650, MDA-MB-231, NMuMG, 3T3-Swiss albino, HCEC, HUVEC, and MCF10A cell lines.

Using the Oris™ Cell Migration Assay offers the following features & benefits:

- **Membrane-free Migration** – perform studies without manipulating transmembrane inserts.
- **Reproducible Results** – obtain well-to-well CV's < 12% due to the unique design.
- **Preserves Cell Morphology** – monitor changes in cell structure in real-time.
- **Versatile** - analyze data using multiple probes in a single well by using a microscope, digital imager, or fluorescence microplate reader.
- **Flexible** - perform kinetic or endpoint cell migration assays without the use of special instrumentation.

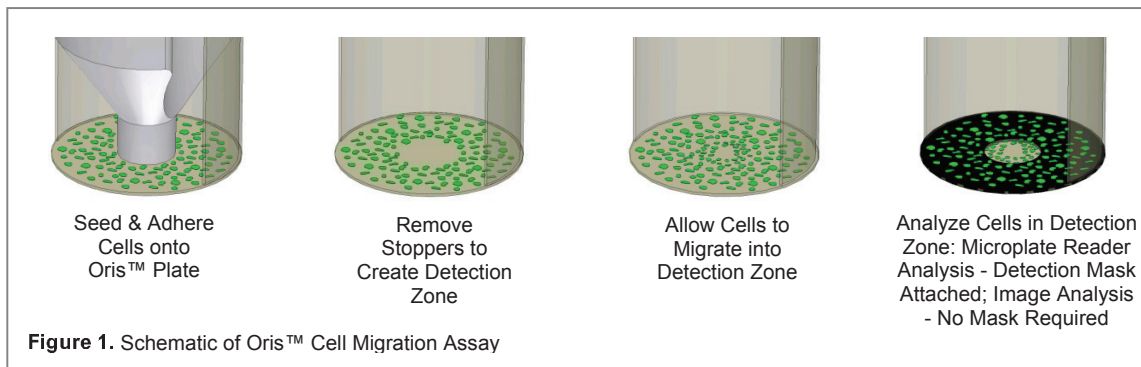

II. ORIS™ PLATE DIMENSIONS

Diameter of Well	6.5 mm
Diameter of Stopper Space (Detection Zone)	2 mm
Suggested Media Volume per Well (populated with Stoppers)	100 µl
Effective Area of Outer Annular Region (seeding region) per Well	30.03 mm ²
Effective Area of Central Detection Zone per Well	3.14 mm ²
Plate Height	14.9 mm
Plate Height with Lid (with Oris™ Cell Seeding Stoppers)	23.9 mm
Offset of Wells (A-1 location, X)	14.4 mm
Offset of Wells (A-1 location, Y)	11.2 mm
Distance between Wells	9 mm (on center)
Well Depth	12.2 mm
Thickness of Well Bottom	0.25 mm
Storage Conditions	Refrigerate (4°C)

Important: Read Instructions Before Performing any Oris™ Assay.

III. MATERIALS PROVIDED

Product No.: CMA1.101

Oris™ 96-well, Tissue Culture Treated (black, clear bottom) Plate with Oris™ Cell Seeding Stoppers, 1
Oris™ Detection Mask, 1
Oris™ Stopper Tool, 1

Product No.: CMA5.101

Oris™ 96-well, Tissue Culture Treated (black, clear bottom) Plates with Oris™ Cell Seeding Stoppers, 5
Oris™ Detection Mask, 1
Oris™ Stopper Tool, 1

IV. MATERIALS REQUIRED

- Biological Cells
- Sterile PBS (containing both Ca⁺⁺ and Mg⁺⁺)
- Complete Cell Culture Growth Medium (containing serum)
- Sterile Pipette Tips/Pipette or Multi-Channel Pipette
- Trypsin or Cell Scraper
- Inverted Microscope (optional)
- Fluorescence Microplate Reader (optional)
- Cell Culture Labeling Medium (phenol red-free/serum-free media)
- Cell Labeling Fluorescent Agent (eg., CellTracker™ Green, Calcein AM)
- *required if performing assay readout via microplate reader.*

Oris™ is a trademark of Platypus Technologies, LLC.

CellTracker™ Green is a trademark of Invitrogen Corporation.

V. CELL MIGRATION ASSAY PROTOCOL

The following steps should be performed in a biological hood using aseptic technique to prevent contamination.

1. Remove the Oris™ Cell Migration Assay Plate from refrigeration and place on lab bench for ~1 hour to allow it to equilibrate to room temperature.
2. Visually inspect the underside of the populated 96-well plate to ensure that the Oris™ Cell Seeding Stoppers are firmly sealed against the bottom of the plate. To inspect the stoppers, turn the plate over and examine the stoppers for sealing (see Figure 2). If incomplete sealing is observed, return the plate to the upright position and use a sterile instrument to gently push the stopper back into the well until sealing is observed.



NOTE: the sealing of the stoppers can be most easily observed if the plate is tipped at an angle and viewed under indirect light to reveal the “bullseye” pattern at the bottom of each well.

3. Apply the Oris™ Detection Mask to the bottom of the 96-well plate if microplate reader data is being collected. The Detection Mask is not necessary if collecting imaging data.

First Time Users: In order to prevent splashing of well contents, familiarize yourself with the attachment and removal of the Detection Mask before any liquids are placed in the wells.

- Orient the chamfered corners of the mask with those of the 96-well plate, ensuring that the A1 corner of the mask is aligned with the A1 well of the plate (see Figure 3).
- Align the holes in the attachment lugs with the bosses on the bottom of the 96-well plate.
- Gently press the mask until it is flush with the bottom of the 96-well plate.



NOTE: It may be necessary to wash the mask with ethanol to remove dust and debris since the mask is **not** sterile. The mask may be applied at any point during the assay. For kinetic assays, it is often most convenient to apply the mask at the beginning of the assay before any liquids are placed in the well. For endpoint assays, using fixed and stained cells, it is often most convenient to apply the mask just before reading assay results.

4. If performing a kinetic analysis of cell migration, pre-label with a fluorescent stain now.
5. Collect cells and prepare a suspension that is 10-fold greater in density than the optimal seeding concentration.

First Time Users: The optimum seeding density of cells must be determined as an integral part of the design of the cell migration assay. Please refer to Appendix I for a discussion of this process.

6. Pipette 100 µl of suspended cells into each test well through one of the side ports of the Oris™ Cell Seeding Stopper.



NOTE: For best results, add or extract media by placing the pipette tip along the wall of the well (see Figure 4). Care should be taken not to disturb the Oris™ Cell Seeding Stopper when introducing the pipette tip into the well. A slender/elongated tip or a gel loading tip may be useful.

7. **IMPORTANT:** Lightly tap the plate on your work surface to evenly distribute well contents (extreme tapping may result in splashing of well contents and lead to contamination).
8. Incubate the seeded plate containing the Oris™ Cell Seeding Stoppers in a humidified chamber (37°C, 5% CO₂) for 4 to 18 hours (cell line dependent) to permit cell attachment.
9. Remove plate from incubator.

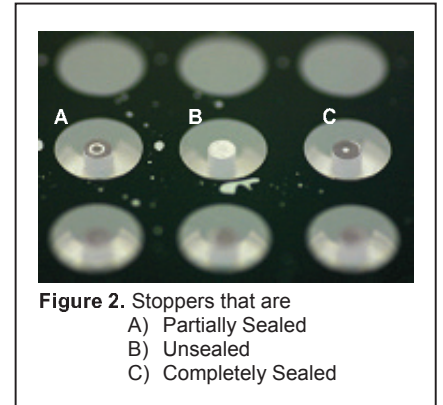


Figure 2. Stoppers that are
A) Partially Sealed
B) Unsealed
C) Completely Sealed

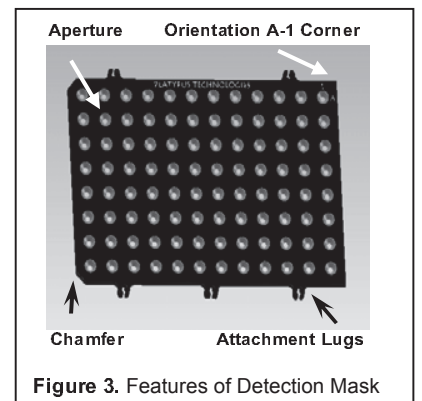


Figure 3. Features of Detection Mask

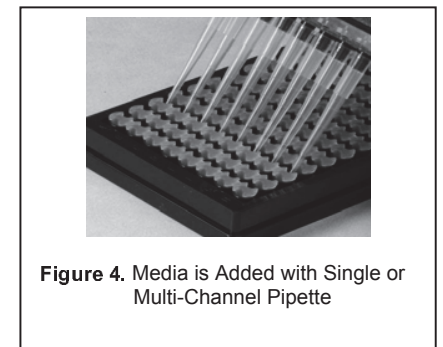


Figure 4. Media is Added with Single or Multi-Channel Pipette

V. CELL MIGRATION ASSAY PROTOCOL, continued

- Designate several 'reference' wells in which the stoppers will remain in place until results are read (t=0 pre-migration controls).
- Using the Oris™ Stopper Tool, remove all other stoppers (see Figure 5).



NOTE: It may be necessary to wash the Oris™ Stopper Tool with 70% ethanol as the Stopper Tool is not sterile.

- Secure the 96-well plate by holding it firmly against the deck of your work space. Slide the tines of the Oris™ Stopper Tool under the backbone of the stopper strip, keeping the underside of the Oris™ Stopper Tool flush with the top surface of the plate.
- Lift the Oris™ Stopper Tool **vertically** to gently remove the stoppers.



NOTE: DO NOT use the Oris™ Stopper Tool as a lever to pry the stoppers from the well (see Figure 5E), as doing so may cause displacement of seeded cells and may distort the detection zone area.

- Remove media with a pipette and **gently** wash wells with 100 µl of sterile PBS (or media) to remove any unattached cells. Do not aspirate using an in-house vacuum.
- Add 100 µl of fresh culture media to each well.
- Incubate plate in a humidified chamber (37°C, 5% CO₂) to permit cell migration. Cells may be examined microscopically throughout the incubation period to monitor progression of migration. Migration time will vary depending upon cell type and experimental design.

- If performing an endpoint analysis of cell migration, stain cells with a fluorescent stain after sufficient migration has occurred. Refer to Section VI and Appendix II for further information on data acquisition and fluorescence staining technique.



NOTE: Oris™ Cell Seeding Stoppers are for single use only; Platypus cannot guarantee the integrity of the stopper material after a second sterilization procedure.

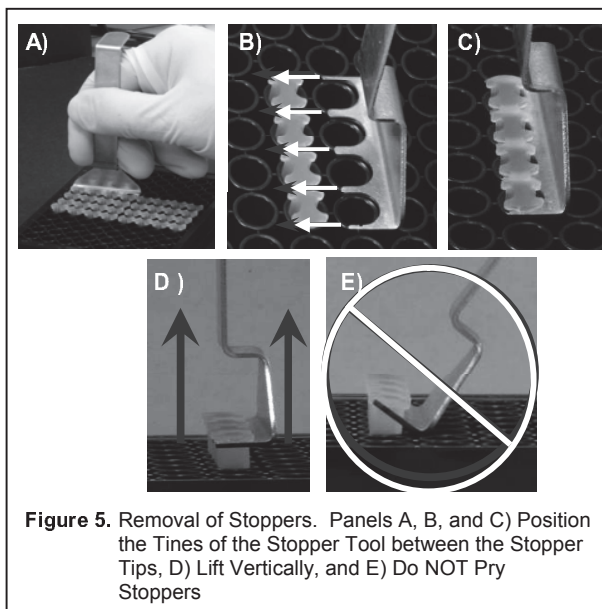


Figure 5. Removal of Stoppers. Panels A, B, and C) Position the Tines of the Stopper Tool between the Stopper Tips, D) Lift Vertically, and E) Do NOT Pry Stoppers

VI. DATA ACQUISITION

The readout of the Oris™ Cell Migration Assay can be conducted at any time, allowing the user to perform a kinetic assay or an endpoint assay. The Oris™ Cell Migration Assay is designed to be used with any commercially available stain or labeling technique. The readout can be performed by using a microscope, a microplate reader, or a High Content Screening or High Content Imaging Analysis platform.

Microscope Analysis

- Cell counting or image capture / analysis software, such as NIH ImageJ freeware, can be used.
- Note: Microscopy observations are possible using phase contrast or bright field microscopy.
- No need to attach the Oris™ Detection Mask to the Oris™ plate.
- Sample data using a colorimetric stain is shown in Figure 6. Wells populated with Oris™ Cell Seeding Stoppers were seeded with 50,000 HT-1080 cells (i.e., 100 μ l of 5×10^5 cells/mL) and incubated overnight. The stoppers were removed from test wells, but remained in place in the pre-migration reference wells until the time of the assay readout. The seeded plate was incubated in a humidified chamber for 24 hours to permit cell migration. Stoppers were removed from the reference wells and all cells were fixed and treated with Wright-Giemsa stain. Images were captured using bright field microscopy and then imported to Image J software for analysis using thresholding. The images below, captured without a detection mask in place, illustrate representative data from pre-migration (t=0 hrs) and post-migration (t=24 hrs) wells (Figure 6A). The graph depicts the average pixel number in the detection zones for each condition (Figure 6B).


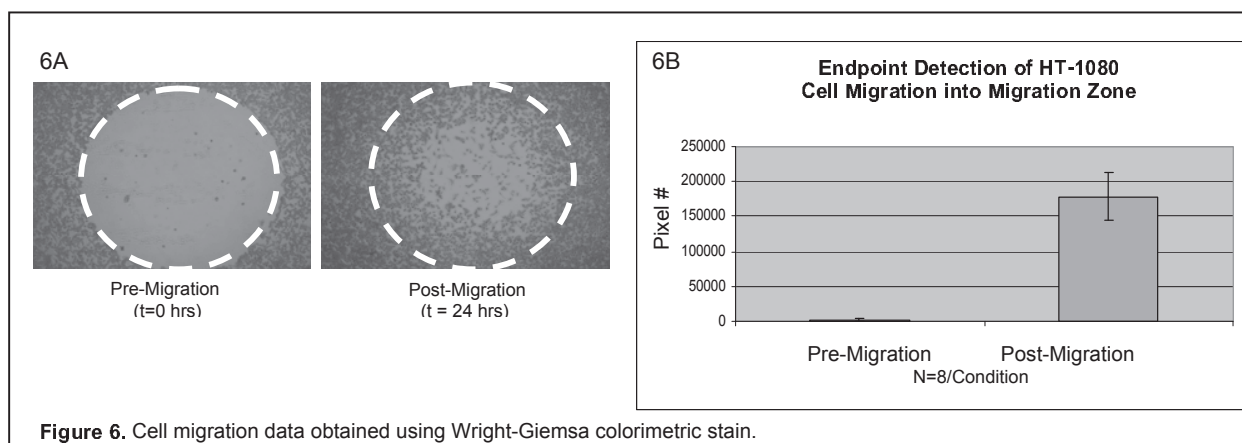


Figure 6. Cell migration data obtained using Wright-Giemsa colorimetric stain.

Microplate Reader Analysis

- Attach the Oris™ Detection Mask to the bottom of the Oris™ plate (refer to Section V, step 3).
- Optimal settings will vary according to the microplate reader make and model. Consult Appendix II and the equipment user manual for your particular instrument.
- The microplate reader MUST be set to read from the bottom of the plate.
- Sample data using a fluorescent stain and microplate reader analysis are shown in Figure 7. Wells populated with Oris™ Cell Seeding Stoppers were seeded with 50,000 HT-1080 cells (i.e., 100 µl of 5×10^5 cells/mL) and incubated overnight. The stoppers were removed from test wells, but remained in place in the pre-migration reference wells until the time of the assay readout. Cells were fluorescently stained with CellTracker™ Green. The seeded plate was incubated in a humidified chamber for 28 hours and at various time points the fluorescence signals in the detection zones were measured using a plate reader. The images below (Figure 7A), captured without a detection mask in place, illustrate representative data from pre-migration (t=0 hrs) and post-migration (t = 21 hrs) wells. The graph depicts a real-time analysis of cell migration that was prepared by transposing the fluorescent signal into cell numbers (Figure 7B).

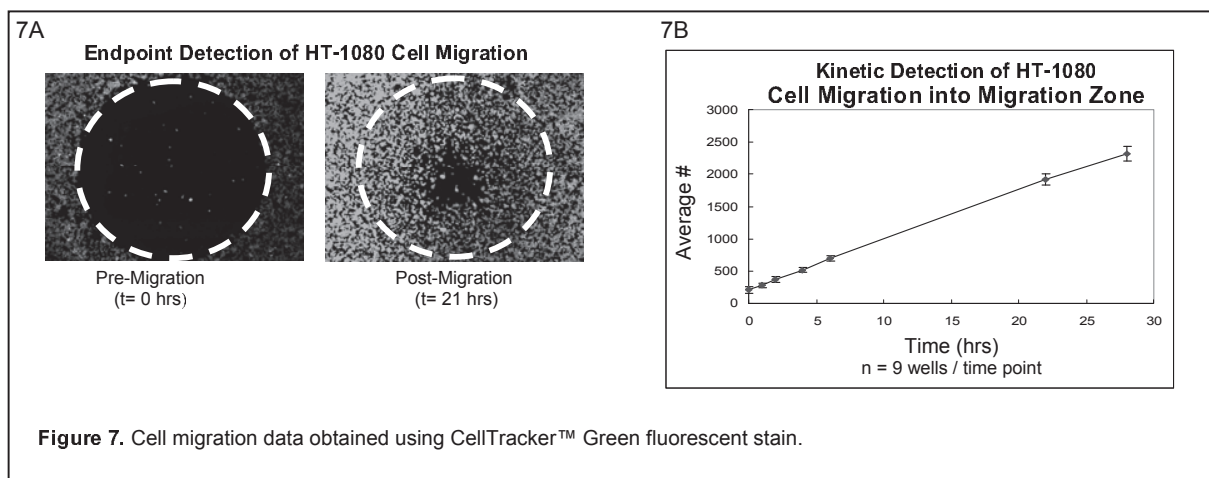

VII. ORDERING INFORMATION

Product Name	Product Number / Product Description	
Oris™ Cell Migration Assays	CMA1.101: 1-pack Oris™ 96-well Plate with Oris™ Cell Seeding Stoppers, 1 Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1	CMA5.101: 5-pack Oris™ 96-well Plates with Oris™ Cell Seeding Stoppers, 5 Oris™ Detection Mask, 1 Oris™ Stopper Tool, 1
Oris™ Cell Migration Assays – Collagen I Coated	CMACC1.101: 1-pack Oris™ 96-well Collagen I Coated Plate with Oris™ Cell Seeding Stoppers, 1 Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1	CMACC5.101: 5-pack Oris™ 96-well Collagen I Coated Plates with Oris™ Cell Seeding Stoppers, 5 Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1
Oris™ Cell Migration Assays – Fibronectin Coated	CMAFN1.101: 1-pack: Oris™ 96-well Fibronectin Coated Plate with Oris™ Cell Seeding Stoppers, 1 Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1	CMAFN5.101: 5-pack Oris™ 96-well Fibronectin Coated Plates with Oris™ Cell Seeding Stoppers, 5 Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1
Oris™ Cell Migration Assays – TriCoated	CMATR1.101: 1-pack Oris™ TriCoated, 96-well Plate with Oris™ Cell Seeding Stoppers, 1 32 wells, Tissue Culture Treated 32 wells, Collagen I coated 32 wells, Fibronectin coated Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1	CMATR5.101: 5-pack Oris™ TriCoated, 96-well Plates with Oris™ Cell Seeding Stoppers, 5 32 wells, Tissue Culture Treated 32 wells, Collagen I coated 32 wells, Fibronectin coated Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1
Oris™ Universal Cell Migration Assembly Kits	CMAU101: 1-pack Oris™-compatible, 96-well Plate, 1 Oris™ Cell Seeding Stoppers, 96 Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1	CMAU505: 5-pack Oris™-compatible, 96-well Plates, 5 Oris™ Cell Seeding Stoppers, 5 x 96 Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1
Oris™ Cell Migration Assembly Kits – Collagen I Coated	CMAUCC1: 1-pack Oris™ Collagen I Coated, 96-well Plate, 1 Oris™ Cell Seeding Stoppers, 96 Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1	CMAUCC5: 5-pack Oris™ Collagen I Coated, 96-well Plates, 5 Oris™ Cell Seeding Stoppers, 5 x 96 Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1
Oris™ Cell Migration Assembly Kit – FLEX	CMAUFL4: 1-kit Oris™-compatible, 96-well Plates, 4 Oris™ Cell Seeding Stoppers, 4 packs of 24* Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1	* sufficient materials for 96 tests
Oris™ Cell Invasion & Detection Assays	CIA101DE: 1-pack Oris™ 96-well plate, 1 Oris™ BME Stock Reagent, 5 mL Oris™ Cell Seeding Stoppers, 96 Oris™ Detection Mask, 1 & Oris™ Stopper Tool, 1 Calcein AM Reagent, 20 µL	CIA200DE: 2-pack Oris™ 96-well plates, 2 Oris™ BME Stock Reagent, 2 x 5 mL Oris™ Cell Seeding Stoppers, 2 x 96 Oris™ Detection Mask, 2 & Oris™ Stopper Tool, 2 Calcein AM Reagent, 2 x 20 µL

VIII. TERMS & CONDITIONS

Certain uses of these products may be covered by U.S. Pat. No. 6,284,197, No. 7,018,838, No. 10,597,118, No. 11,342,413, and No. 60,836,109, licensed to PLATYPUS. Certain applications of PLATYPUS products may require licenses from other parties. Determining the existence and scope of such third party intellectual property is the responsibility of the PURCHASER. Purchase of the product provides the PURCHASER with a limited non-transferable license under any PLATYPUS patents or patent applications to use the product for internal research unless there is a written limitation to this license in the product literature. PURCHASER is responsible for carefully reviewing the product literature and respecting any limitations to this license, e.g. limitations for commercial use or research by for-profit institutions. These products may not be resold, modified for resale, used to manufacture commercial products, or used to develop commercial products without the express written approval of PLATYPUS.

These products are intended for research or laboratory use only and are not to be used for any other purposes, including, but not limited to, unauthorized commercial purposes, *in vitro* diagnostic purposes, *ex vivo* or *in vivo* therapeutic purposes, investigational use, in foods, drugs, devices or cosmetics of any kind, or for consumption by or use in connection with or administration or application to humans or animals. PLATYPUS warrants that its products shall conform substantially to the description of such goods as provided in product catalogues and literature accompanying the goods until their respective expiration dates or, if no expiration date is provided, for 6 months from the date of receipt of such goods. PLATYPUS will replace, free of charge, any product that does not conform to the specifications. This warranty limits PLATYPUS's liability only to the replacement of the nonconforming product.

THIS WARRANTY IS EXCLUSIVE AND PLATYPUS MAKES NO OTHER WARRANTY, EXPRESS OR IMPLIED, INCLUDING WITHOUT LIMITATION, ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. The stated express warranties, and the remedy provided for breach thereof, are in lieu of all other liability or obligations of PLATYPUS for any damages whatsoever arising out of or in connection with the delivery, use, misuse, performance, or the inability to use any of its products. IN NO EVENT SHALL PLATYPUS BE LIABLE UNDER ANY LEGAL THEORY (INCLUDING BUT NOT LIMITED TO CONTRACT, NEGLIGENCE, STRICT LIABILITY IN TORT, OR WARRANTY OF ANY KIND) FOR ANY INDIRECT, SPECIAL, INCIDENTAL, CONSEQUENTIAL, OR EXEMPLARY DAMAGES (INCLUDING BUT NOT LIMITED TO LOST PROFITS) EVEN IF PLATYPUS HAD NOTICE OF THE POSSIBILITY OF SUCH DAMAGES. Without limiting the effect of the preceding sentence, PLATYPUS's maximum liability, if any, shall not exceed the purchase price paid by PURCHASER for the product.

This warranty shall not be effective if PLATYPUS determines, in its sole discretion that PURCHASER has altered or misused the goods or has failed to use or store them in accordance with instructions furnished by PLATYPUS. PLATYPUS's sole and exclusive liability and PURCHASER's exclusive remedy with respect to goods proved to PLATYPUS's satisfaction (applying analytical methods reasonably selected by PLATYPUS) to be defective or nonconforming shall be the replacement of such goods free of charge, upon the return of such goods in accordance with our instructions, although at its discretion, PLATYPUS may provide a credit or refund. If PLATYPUS manufactures custom goods for PURCHASER based on instructions, specifications, or other directions provided by PURCHASER, PLATYPUS shall not be liable for the lack of sufficiency, fitness for purpose or quality of the goods to the extent attributable to such instructions, specifications, or other directions. PLATYPUS shall not be liable for any loss, damage or penalty as a result of any delay in or failure to manufacture, deliver or otherwise perform hereunder due to any cause beyond PLATYPUS's reasonable control.

PLATYPUS shall not be liable for injury or damages resulting from the use or misuse of any of its products.

APPENDIX I: Determining Optimal Cell Seeding Concentration

This procedure is intended to assist in determining the cell seeding density needed to achieve confluency of your cell line when using the Oris™ Cell Migration Assay. The intended goal is to achieve 90-95% confluency of the monolayer surrounding the Oris™ Cell Seeding Stoppers without overgrowth.

1. A suggested starting point is to evaluate three serial dilutions at the cell densities shown below. The cell seeding area of the well with the stopper in place is ~ 0.3 cm². Based on the typical seeding density of your particular cell line, you can infer a different cell number for your first serial dilution and adjust the numbers below accordingly.
2. Prepare a log-phase culture of the cell line to be tested. Collect cells and determine the total number of cells present.
3. Pellet cells by centrifugation. Prepare three serial dilutions at final concentrations of 1.0 x 10⁶, 0.5 x 10⁶ and 0.25 x 10⁶ cells/mL.
4. Dispense 100 µl of cell suspension per well into the 96-well plate to result in the following plate layout:

Column	1	2	3
Cells / well	100,000	50,000	25,000
Number of wells	8	8	8

5. Incubate the plate in a humidified chamber (37°C, 5% CO₂) for 4 - 18 hours (cell line dependent) with cell seeding stoppers in place to allow the cells to firmly attach to the well surface.
6. Following cell attachment, remove the Oris™ Cell Seeding Stoppers from each well (see Figure 5) and **gently** wash the wells with PBS to remove non-attached cells.
 - Secure the 96-well plate by holding it firmly against the deck of your work space. Slide the tines of the Oris™ Stopper Tool under the backbone of the stopper strip, keeping the underside of the tool flush with the top surface of the plate.
 - Lift the Oris™ Stopper Tool **vertically** to gently remove the stopper. Do not use the Oris™ Stopper Tool as a lever to pry the stoppers from the well as doing so may cause displacement of the seeded cells.
7. Without a Detection Mask in place, use a microscope to visually inspect each well to determine the minimum cell seeding concentration that yields a confluent monolayer at the perimeter of the detection zone.

At this point, if you plan to obtain the results of the Oris™ Cell Migration Assay via colorimetric or microscopy analysis, you have successfully determined the optimal cell seeding concentration to be used in Step 5 of the Cell Migration Assay Protocol.

APPENDIX II: Determining Optimal Fluorescence Microplate Reader Settings

This procedure is intended to assist in optimizing your instrument settings when using a fluorescence microplate reader to capture data from the Oris™ Cell Migration Assay.

1. Using the optimal cell seeding concentration determined in Appendix I, perform a cell migration assay per Section V, Cell Migration Assay Protocol using culture conditions expected to result in robust cell migration. Be sure to include equal numbers of pre-migration reference wells (stoppers left in place until staining) and post-migration test wells (stoppers removed after cell attachment period). A minimum of 8 wells per condition are recommended.
2. Perform the desired fluorescent staining technique.

The Oris™ Cell Migration Assay has been designed to work with all types of fluorescent stains and staining techniques. The precise method for staining cells with fluorescent stains varies according to the nature of the individual stain. It is important to stain cells using a fluorescent reagent that uniformly stains cells. Probes affected by experimental conditions will increase variability of results and reduce correlation between fluorescence signal and cell migration. Please consult the manufacturer of your fluorescent stain for specific considerations.

The following is an example Fluorescent Staining Protocol for using Calcein AM:

- a. To stain one fully-seeded 96-well plate, combine 5 µl of Calcein AM (1 mg/mL in dry DMSO) with 10 mL of phenol red-free and serum-free media or 1x PBS (containing both Ca⁺⁺ and Mg⁺⁺). Protect diluted Calcein AM solution from light until ready to use in step d.
 - b. Carefully remove culture medium from wells.
 - c. Wash wells with 100 µl of PBS (containing both Ca⁺⁺ and Mg⁺⁺).
 - d. Add 100 µl of diluted Calcein AM solution to each well.
 - e. Incubate plate at 37°C for 30 - 60 minutes.
 - f. Attach mask and read promptly with microplate reader using appropriate filter set and sensitivity/gain settings (for a BioTek Synergy™ HT microplate reader, use 485/528 nm excitation/emission filters, sensitivity 55 nm).
3. If not already in place, apply the Oris™ Detection Mask to the plate. Using the bottom probe of a fluorescence microplate reader, obtain the fluorescence reading from each well. To achieve the optimal dynamic range, adjust the instrument settings (e.g., gain) to result in the greatest difference in fluorescence signal between pre-migration and post-migration wells. Refer to the instrument manual for your microplate reader for further guidance on instrument settings.

You have now successfully determined the optimal cell seeding concentration (to be used in Step 5 of the Cell Migration Assay Protocol) and microplate reader settings for analysis of cell migration using a fluorescence microplate reader.

AMSBIO | www.amsbio.com | info@amsbio.com



UK & Rest of the World
184 Park Drive, Milton Park
Abingdon OX14 4SE, U.K.
T: +44 (0) 1235 828 200
F: +44 (0) 1235 820 482



North America
1035 Cambridge Street,
Cambridge, MA 02141.
T: +1 (617) 945-5033 or
T: +1 (800) 987-0985
F: +1 (617) 945-8218



Germany
Bockenheimer Landstr. 17/19
60325 Frankfurt/Main
T: +49 (0) 69 779099
F: +49 (0) 69 13376880



Switzerland
Via Lisano 3
(CP.683)
6900 Massagno
T: +41 (0) 91 604 55 22
F: +41 (0) 91 605 17 85