BioVision rev. 09/17 For research use only

Active Human Cathepsin L

CATALOG #: 1135-5 5 μg 1135-100 100 μα

1135-100 100 μg 1135-1000 1 mg

SOURCE: E. coli

ALTERNATIVE NAMES: CTSL, CTSL1

FORM: Liquid.

FORMULATION: 0.5 mg/ml in 25 mM Sodium Acetate, 50 mM NaCl, pH 5.5

containing 50% glycerol

PURITY: ≥90% by SDS-PAGE

MOL. WT.: 24.2 kDa (114-333 aa)

STORAGE CONDITIONS: Store at -80° C, stable for at least 6 months as supplied. Use

as provided or dilute with 50% glycerol to 0.1 mg/ml. Aliquot immediately and store at -80°C. Avoid freeze/thaw cycles.

APPLICATIONS: Recombinant Cathepsin L is useful in inhibitor screening,

kinetic studies and cleavage of substrates in vitro or in tissue

culture.

ACTIVITY ASSAY: BioVision's Active Cathepsin L has been activated through

auto-catalytic activation of Procathepsin L (Cat # 7581) at low pH. The enzyme is fully active as seen from its ability to cleave a fluorogenic substrate Z-Phe-Arg-AFC (K142-100).

SPECIFIC ACTIVITY: > 2,000 mU/mg (1 U = digestion of 1 μmole of Z-Phe-Arg-AFC

substrate (K142-100) per minute at 37°C in 50 mM Sodium

Acetate, 100 mM NaCl, pH 5.5).

BACKGROUND: Cathepsin L is a lysosomal protease that can be secreted from


the cell and is implicated to play a role in protein degradation,

arthritis, apoptosis, and cancer.

REFERENCES:


- Joseph, L.J., et al. (1988) J. Clin. Invest. 81:1621-1626.
- Turk, V., et al. (2001) EMBO J. 20:4629-4634.

For Research Use Only! Not to be used in humans.

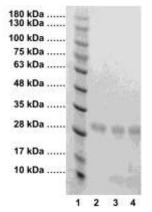

SDS-PAGE (4-20%) of h-Cathepsin L:

- 1: Protein Marker
- 2: Human Cathepsin L (5 µg)
- 3: Human Cathepsin L (8 µg)
- 4: Human Cathepsin L (10 µg)

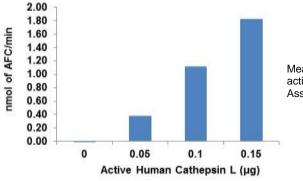

Measurement of h-Cathepsin L activity by Cathepsin L Activity Assay Kit (K142-100)

RELATED PRODUCTS:

- Human CellExp™ Cathepsin B, human recombinant (Cat. No. 7408-10)
- Human CellExp™ Cathepsin D, human recombinant (Cat. No. 7409-10)
- Human CellExp™ Cathepsin S, human recombinant (Cat. No. 7277-10, -50, -1000)
- Human CellExp™ Cathepsin L1, human recombinant (Cat. No. 7410-10)
- Cathepsin K, Active, human recombinant (Cat. No. 7600-5, -50)
- Cathepsin K, Active, mouse recombinant (Cat. No. 7597-5, -50)
- Cathepsin K, Active, rat recombinant (Cat. No. 7598-5, -50)
- Cathepsin B (active, human) (Cat. No. 1021-5)
- Cathepsin D (active, human) (Cat. No. 1022-5)
- Cathepsin H (active, human) (Cat. No. 1023-5)
- Procathepsin K, human recombinant (Cat. No. 1026-10, -50, -1000)
- Procathepsin K, mouse recombinant (Cat. No. 1027-10, -50, -1000)
- Procathepsin K, rat recombinant (Cat. No. 1029-10, -50, -1000)