

Catalog Number: IR13-36

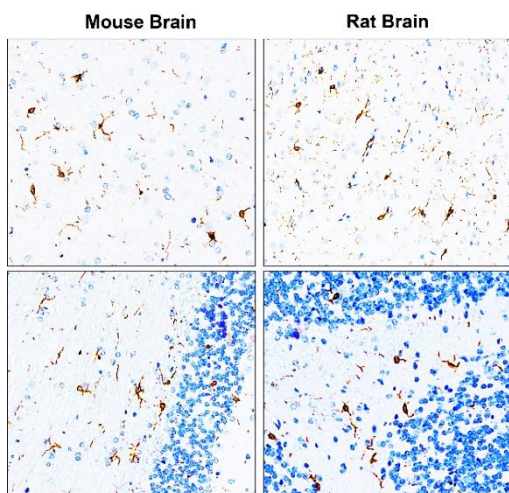
Anti-Iba1 / AIF antibody

Package size: 25, 100 µl
Store at: -20°C
MW (kDa): 15

Overview

Product Name	Iba1 / AIF antibody
Product Number	IR13-36
Gene Description	allograft inflammatory factor 1
Clonality	Polyclonal
Host	Rabbit
Species Reactivity	Human, mouse, rat
Recommended Applications Dilutions	Western Blot 1:500 – 1:1000 Immunofluorescence 1:100 – 1:300 Immunohistochemistry (Paraffin) 1:100 – 1:300 Immunohistochemistry (Frozen) 1:100 – 1:300
Storage Buffer	100mM Tris Glycine, 1% BSA, 20% Glycerol (pH7). 0.025% ProClin 300 was added as a preservative
Concentration	0.7 mg/ml
Purity	Affinity column purified
Storage	Store at +4°C for short term storage. Long time storage is recommended at -20°C
Notes	Gently mix before use. Optimal concentrations and conditions for each application should be determined by the user.

IR13-36 anti-Iba1/AIF1 antibody IHC image

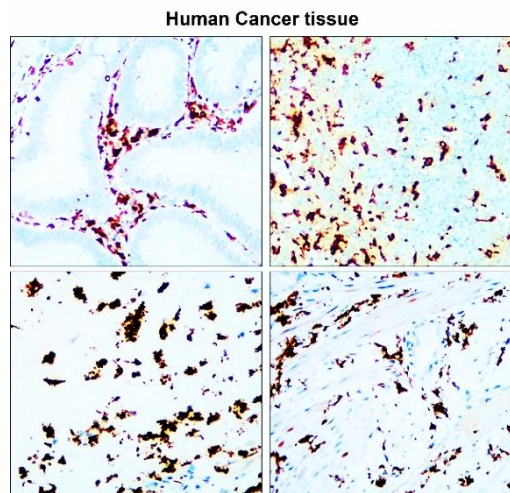


Immunohistochemical analysis of paraffin embedded Mouse / Rat brain tissue labeling Iba1 with IR13-36 at 1/200.

Epitope Retrieval methods: Citrate Buffer, pH 6.0

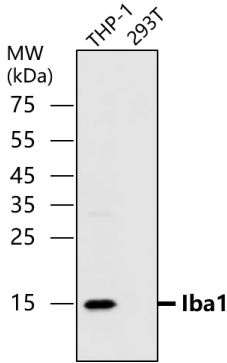
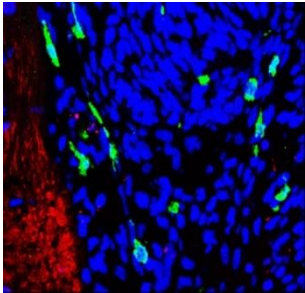
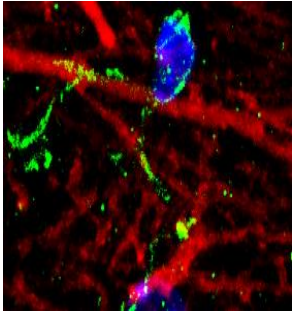
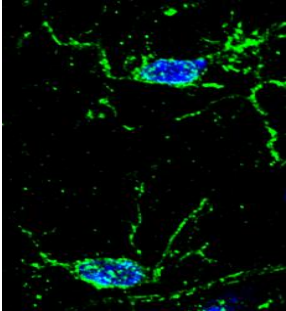
Data from: [Li-Tzung Biotechnology. Co](#)

IR13-36 anti-Iba1/AIF1 antibody IHC image



Immunohistochemical analysis of paraffin embedded Human cancer tissue labeling Iba1 with IR13-36 at 1/200.

Epitope Retrieval methods: Citrate Buffer, pH 6.0

IR13-36 anti-Iba1/AIF1 antibody WB image	Customer feedback Image IR13-36 anti-Iba1/AIF1 antibody IHC-Fr image
Iba1 Microglia marker 	
Iba1/AIF1 antibody - at 1/500 dilution Lysates at 45 µg per lane, This blot was produced using a 15% SDS-PAGE. PVDF membrane was then blocked with 3%BSA for an hour before being incubated with IR13-36 overnight at 4°C.	Frozen sectioned for immunohistochemical analysis. Sample: Frozen sectioned Rat brain. Green: Iba1 antibody (IR13-36) Red: TUJ1 antibody Blue: DAPI
Customer feedback Image IR13-36 anti-Iba1/AIF1 antibody IHC-Fr image	Customer feedback Image IR13-36 anti-Iba1/AIF1 antibody IHC-Fr image
	
Immunohistochemical of frozen sections. Sample : mouse cerebellum. Green : Iba1 antibody (IR13-36) : 1-200 Anti-rabbit 488 : 1-500. DAPI (blue) was used as the nuclear counter stain.	Immunohistochemical of frozen sections. Sample : mouse hippocampus. Green : Iba1 antibody (IR13-36) : 1-200 Anti-rabbit 488 : 1-500. DAPI (blue) was used as the nuclear counter stain.

IReal Biotechnology Co., Ltd.



**ISO 13485:2016 Quality
Management System
Certified**

