


## Minocycline hydrochloride

<b>Cat. No.:</b>	HY-17412
<b>CAS No.:</b>	13614-98-7
<b>Molecular Formula:</b>	C <sub>23</sub> H <sub>28</sub> ClN <sub>3</sub> O <sub>7</sub>
<b>Molecular Weight:</b>	493.94
<b>Target:</b>	Bacterial; Antibiotic; HIF/HIF Prolyl-Hydroxylase; Apoptosis; MDM-2/p53; Potassium Channel; Calcium Channel
<b>Pathway:</b>	Anti-infection; Metabolic Enzyme/Protease; Apoptosis; Membrane Transporter/Ion Channel; Neuronal Signaling
<b>Storage:</b>	4°C, sealed storage, away from moisture * In solvent : -80°C, 1 years; -20°C, 6 months (sealed storage, away from moisture)

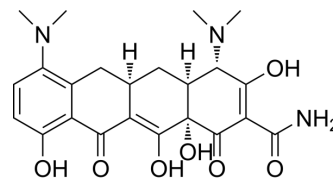

HCl

### SOLVENT & SOLUBILITY

<b>In Vitro</b>	DMSO : 19.23 mg/mL (38.93 mM; Need ultrasonic)			
	H <sub>2</sub> O : 9.09 mg/mL (18.40 mM; Need ultrasonic)			
		Solvent Concentration	Mass	
			1 mg	5 mg
			10 mg	
<b>Preparing Stock Solutions</b>	1 mM	2.0245 mL	10.1227 mL	20.2454 mL
	5 mM	0.4049 mL	2.0245 mL	4.0491 mL
	10 mM	0.2025 mL	1.0123 mL	2.0245 mL
Please refer to the solubility information to select the appropriate solvent.				
<b>In Vivo</b>	1. Add each solvent one by one: PBS Solubility: 7.69 mg/mL (15.57 mM); Clear solution; Need ultrasonic			

### BIOLOGICAL ACTIVITY

<b>Description</b>	Minocycline hydrochloride is an orally active, potent and BBB-penetrated semi-synthetic tetracycline antibiotic. Minocycline hydrochloride is a hypoxia-inducible factor (HIF)-1 $\alpha$ inhibitor. Minocycline hydrochloride shows anti-cancer, anti-inflammatory, and glutamate antagonist effects. Minocycline hydrochloride reduces glutamate neurotransmission and shows neuroprotective properties and antidepressant effects. Minocycline hydrochloride inhibits bacterial protein synthesis through binding with the 30S subunit of the bacterial ribosome, resulting in a bacteriostatic effect <sup>[1][2][3][4][5][6][7]</sup> .
<b>IC<sub>50</sub> &amp; Target</b>	L-type calcium channel
<b>In Vitro</b>	Minocycline hydrochloride (0-100 $\mu$ M, 24-72 h) suppresses proliferation and clonogenic activity of ovarian cancer cell-lines (OVCAR-3, SKOV-3 and A2780) <sup>[3]</sup> . Minocycline hydrochloride (0-100 $\mu$ M, 24-48 h) arrests cell cycle through inhibition of cyclins and suppression of DNA

incorporation<sup>[3]</sup>.

Minocycline hydrochloride (0-100  $\mu$ M, 72 h) induces cell apoptosis in ovarian cancer cell lines<sup>[3]</sup>.

Minocycline hydrochloride shows direct neuronal protection, and this mode of protection is likely to be associated with the preservation of mitochondrial integrity and cytochrome c, followed by the suppression of caspase-dependent as well as caspase-independent cell death<sup>[2]</sup>.

Minocycline hydrochloride leads to suppression of Hypoxia-inducible factor (HIF)-1 $\alpha$  accompanied by up-regulation of p53 protein levels and inactivation of AKT/mTOR/p70S6K/4E-BP1 pathway<sup>[6]</sup>.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

#### Cell Proliferation Assay<sup>[3]</sup>

Cell Line:	Human ovarian cancer cell lines (OVCAR-3, SKOV-3 and A2780) and primary cells (HEK-293, HMEC, HUVEC, ATCC)
Concentration:	0, 1, 10, 50 and 100 $\mu$ M
Incubation Time:	24, 48 or 72 h
Result:	Inhibited proliferation of OVCAR-3, SKOV-3 and A2780 cells in a concentration-dependent manner, with IC <sub>50</sub> values of 62.0, 56.1 and 59.5 $\mu$ M, respectively. Had no effect on the viability of HEK-293 or HUVEC.

#### Western Blot Analysis<sup>[3]</sup>

Cell Line:	OVCAR-3, SKOV-3 and A2780 cells
Concentration:	0, 10, 50 and 100 $\mu$ M
Incubation Time:	72 h
Result:	Expressed lower levels of cyclins A, B and E. Increased caspase-3 levels by more than 3.0 fold in the 100 $\mu$ M. Minocycline-activated caspase-3 in turn led to cleavage of PARP-1. Increased the degradation product p89 of PARP-1 by caspase-3.

#### Cell Cycle Analysis<sup>[3]</sup>

Cell Line:	OVCAR-3, SKOV-3 and A2780 cells
Concentration:	0, 10, 50 and 100 $\mu$ M
Incubation Time:	24 or 48 h
Result:	Arrested cells in the G0-G1 phase in a concentration and time-dependent manner. Declined percentage of cells in the S and G2-M phases in excess of 80% each at 100 $\mu$ M.

#### In Vivo

Minocycline hydrochloride (0-30 mg/kg, orally, daily for 4 weeks) suppresses OVCAR-3 tumor growth in female nude mice<sup>[3]</sup>.

Minocycline hydrochloride (IP) is an effective neuroprotective agent in animal models of cerebral ischemia when given in high doses intraperitoneally<sup>[1]</sup>.

Minocycline hydrochloride (0-40 mg/kg, IP, once) significantly attenuates METH-induced hyperlocomotion and the development of behavioral sensitization in mice<sup>[2]</sup>.

Minocycline hydrochloride (3 and 10 mg/kg, IV, once) is effective at reducing infarct size in a Temporary Middle Cerebral Artery Occlusion model (TMCAO)<sup>[1]</sup>.

Minocycline hydrochloride (3-10 mg/kg, IV, once) results in serum levels (at 3 mg/kg) similar to that achieved in humans after a standard 200 mg dose<sup>[1]</sup>.

Minocycline hydrochloride attenuates ischemia-induced ventricular arrhythmias in rats. This effect may be associated with activations of PI3K/Akt signaling pathway, mitochondrial KATP channels and L-type Ca<sup>2+</sup> channels<sup>[7]</sup>.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

Animal Model:	Female nude mice (6 weeks old, 9 per group, OVCAR-3 cells were injected s.c. into the left flank of each mouse) <sup>[3]</sup>
Dosage:	10 or 30 mg/kg
Administration:	Administered orally in the drinking water, initiated on day 8 of cell inoculation, daily for 4 weeks
Result:	Suppressed OVCAR-3 tumor growth in these female nude mice, and reduced microvessel density.
Animal Model:	Male Balb/cAnNCrIcrIj mice (8 weeks old, 23-30 g, methamphetamine (METH, 3 mg/kg) was injected subcutaneously (s.c.) in a volume of 10 ml/kg) <sup>[2]</sup>
Dosage:	0, 10, 20, or 40 mg/kg
Administration:	IP, once, 30 min before the administration of METH
Result:	Significantly attenuated METH-induced hyperlocomotion and the development of behavioral sensitization in mice at 40 mg/kg. Did not exert any effect on the induction of METH-induced hyperthermia in mice. Significantly attenuated the reduction of DA and DOPAC in the striatum. Significantly attenuated the reduction of DAT-immunoreactivity in the mouse striatum. Significantly attenuated the increase in MAC1-immunoreactivity in the striatum after the administration of METH.
Animal Model:	Male Sprague-Dawley rats (270-330 g, TMCAO model) <sup>[1]</sup>
Dosage:	3 mg/kg and 10 mg/kg
Administration:	IV, once, 4, 5, or 6 hours post TMCAO
Result:	Reduced infarct size by 42% while 10 mg/kg reduced infarct size by 56% at doses of 3 mg/kg; significantly reduced infarct size at 5 hours by 40% at doses of 10 mg/kg and the 3 mg/kg dose significantly reduced infarct size by 34%. With a 6 hour time window there was a non-significant trend in infarct reduction.
Animal Model:	Male Sprague-Dawley rats (270-330 g) <sup>[1]</sup>
Dosage:	3, 10, or 20 mg/kg
Administration:	IV, once
Result:	Peak concentrations of serum levels of minocycline averaged 3.6, 13.0 and 28.8 mg/L with 3, 10 and 20 mg/kg doses respectively. The serum levels of minocycline at a 3 mg/kg dose (3.6 mg/L) were similar to that reported in humans after a standard 200 mg dose. Did not significantly affect hemodynamic and physiological variables.

## CUSTOMER VALIDATION

- Adv Funct Mater. 2023 Jul 7.

- Nat Commun. 2022 Mar 2;13(1):1116.
- Brain Behav Immun. 2021 Jan;91:519-530.
- Brain Behav Immun. 2021 Jan;91:505-518.
- Brain Behav Immun. 2020 Oct;89:400-413.

See more customer validations on [www.MedChemExpress.com](http://www.MedChemExpress.com)

## REFERENCES

---

- [1]. Xu L, et al. Low dose intravenous minocycline is neuroprotective after middle cerebral artery occlusion-reperfusion in rats. BMC Neurol. 2004 Apr 26;4:7.
- [2]. Zhang L, et al. Protective effects of minocycline on behavioral changes and neurotoxicity in mice after administration of methamphetamine. Prog Neuropsychopharmacol Biol Psychiatry. 2006 Dec 30;30(8):1381-93.
- [3]. Pourgholami MH, et al. Minocycline inhibits growth of epithelial ovarian cancer. Gynecol Oncol. 2012 May;125(2):433-40.
- [4]. Molina-Hernández M, et al. Antidepressant-like actions of minocycline combined with several glutamate antagonists. Prog Neuropsychopharmacol Biol Psychiatry. 2008 Feb 15;32(2):380-6.
- [5]. Ritchie DJ, et al. A review of intravenous minocycline for treatment of multidrug-resistant Acinetobacter infections. Clin Infect Dis. 2014 Dec 1;59 Suppl 6:S374-80.
- [6]. Ataie-Kachoei P, et al. Minocycline attenuates hypoxia-inducible factor-1 $\alpha$  expression correlated with modulation of p53 and AKT/mTOR/p70S6K/4E-BP1 pathway in ovarian cancer: in vitro and in vivo studies. Am J Cancer Res. 2015 Jan 15;5(2):575-88.
- [7]. Hu X, Wu B, Wang X, Xu C, He B, Cui B, Lu Z, Jiang H. Minocycline attenuates ischemia-induced ventricular arrhythmias in rats. Eur J Pharmacol. 2011 Mar 11;654(3):274-9.
- 

**Caution: Product has not been fully validated for medical applications. For research use only.**

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: [tech@MedChemExpress.com](mailto:tech@MedChemExpress.com)

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA